O₂MODA

Eye Care

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Eye Beauty Massager

Daily eye-area care

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— 、 Safety Warnings

1 、 During use, you may feel a kneading or massaging sensation around the eyes and eyelids. There may also be slight pressure in the nasal cavity or upper jaw. There should be no sharp or stinging sensation.

If pain occurs, it may be due to damaged skin contact—stop immediately.

If caused by high intensity, reduce the output level.

2 、 Sensation may differ between left and right eyes. This does not affect performance and may be due to:

- a. Different electrode contact tightness
- b. Different skin moisture levels
- c. Natural biological differences between left and right eyes

3 、 Read the manual carefully before use.

4 、 Electrodes must be used only in designated areas.

5 、 Do not drop, expose to heat, water, or corrosive environments.

6 、 Hold the plug when connecting/disconnecting charging cable.

7 、 Do not use near microwave devices.

8 、 If you have a pacemaker or electronic implant, consult a doctor before use.

9 、 Do not disassemble or modify the device.

10 、 Not intended for use by individuals (including children) with reduced physical, sensory, or mental capabilities unless supervised.

11、Children must be supervised to prevent misuse.

12、Use only with 5V  low voltage safety power supply.

二、Disclaimer

The company makes no warranties, including fitness for a particular purpose.

Not responsible for errors or damages arising from use.

This manual contains proprietary information protected by law.

No reproduction or translation without permission.

Content may change without notice.

三、Overview

This product is an EMS Eye Massager that uses microcurrent technology for beauty care.

Uses low intensity microcurrent to simulate massage

Helps relax the skin around the eyes

Gently revitalizes and reduces fatigue

Provides a comfortable and soothing experience

The microcurrent is specially programmed to ensure a smooth and gentle experience without noticeable tingling.

Conductive silicone pads deliver current to the eye area, simulating massage and restoring a refreshed feeling.

四、Functions & Features

1、Functions

Suitable for daily eye care

Helps refresh and relax eye area after prolonged use

Provides comfort and cosmetic care experience

2、Features

- **EMS Microcurrent Technology**

Simulates gentle finger massage for daily care

- **Ergonomic Design**

Conductive silicone fits evenly around the eye area

(Does not contact the eyeball)

- **Revitalizing Effect**

Gentle stimulation helps restore vitality

- **Multiple Modes**

Pulse mode & Continuous mode

- **Modern Integrated Design**

Stylish appearance with interchangeable lenses

五、Contraindications

Do NOT use if you have:

Damaged skin around eyes

Acute eye inflammation

Conjunctivitis / keratitis / corneal ulcers

Glaucoma

Epilepsy

Brain vascular stents

Head tumors

Pacemaker users

Bleeding disorders

Severe heart disease

Allergic constitution

Pregnancy

This product is for external cosmetic massage only.

Consult a doctor for medical conditions.

六、Technical Specifications

(This product features 20 intensity levels that adjust only the output current.

The pulse frequency, pulse width, and waveform remain fixed and unchanged.)

Power Supply: DC 3.7V (built-in lithium battery)

Pulse Frequency: 20Hz $\pm 10\%$

Pulse Width: 300 μs $\pm 10\%$

Waveform: Square wave


Output Current: $\leq 1.5\text{mA}$

Output Power: 0.2A

Carrier Frequency: 3kHz $\pm 10\%$

Carrier Pulse Width: 100 μs $\pm 10\%$

七、Product Description

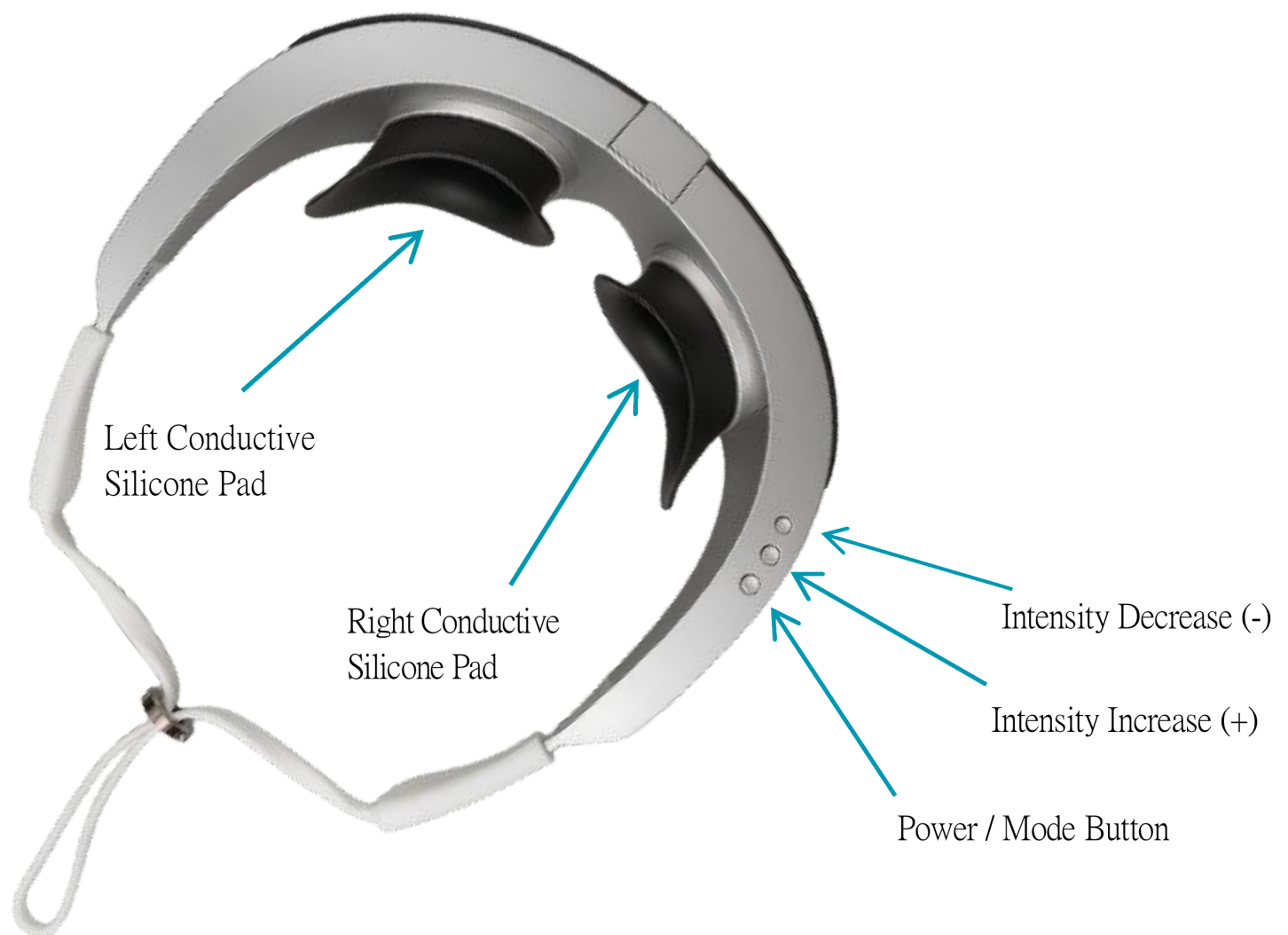八、Installation & Usage

1、Charging

- Connect adapter to device and power outlet
- Red light = charging
- Light off = fully charged
- Flashing:
 - 2 times/sec = normal operation
 - 5 times/sec = low battery

2. Wearing the Device

- Moisten skin around eyes before use
- Attach silicone pads firmly
- Adjust strap for comfortable fit
- Avoid excessive moisture or sweat

3. Operation

① **Power On** : Long press → voice prompt: “Power On”

*Default: Pulse Mode, Level 1

② **Mode Switch** : Short press to switch:

- Pulse Mode
- Continuous Mode

③ **Intensity Adjustment** : 20 levels total

- Levels 1 – 10: low
- Levels 11 – 20: high

Beginners should start at low levels and increase gradually.

④ **Timer** : Default: 15 minutes

Auto voice: “Session complete”

Auto shutdown

⑤ **Pause Function** : Press both +/- buttons → pause

Press any button → resume gradually

⑥ **Power Off** : Long press → “Power Off”

九、Warranty

- 1-year free warranty
- Warranty does NOT cover:
 - a.Unauthorized modifications
 - b.Damage from drops, water, misuse
 - c.Incorrect operation
- Proof of purchase required
- Out-of-warranty repairs may incur fees
- Service via authorized centers
- Consumables & accessories not covered

+ 、Package Contents

- Main Unit ×1
- USB Charging Cable ×1
- Elastic Strap ×1
- User Manual ×1
- Warranty Card ×1

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